

WIRELESS HANDHELD ULTRASOUND SCANNERS



**Available in both Black & White and Color Doppler
with Linear, Convex, Micro Convex, and Transvaginal Probes**

The **High Frequency Linear Probes** work well for close to the surface scans such as tendons on horses and small animal pregnancy with penetration depths of up to 100 mm. The **Lower Frequency Convex & Micro Convex Probes** work best for abdominal organs such as kidneys, liver, spleen, and heart. They also work for larger animal pregnancy such as sheep, goats, lamas, and alpacas with penetration depths of up to 305 mm.

Features:

- 2 Year Warranty
- Built-in rechargeable battery with 3 to 5 hours of battery life
- Battery charge time – 2 hours
- Comes with wireless battery charger
- Images stored in Photo Galley can be attached to emails and texts
- Image Adjustments – BGain, TGC, DYN, Focus, Depth, Harmonic, Denoise, Color Gain, Steer, & PRF
- Working system: Apple iOS and Windows
- Measure: Length, Area, Angle, Heart Rate, Obstetrics
- Cineplay: auto or manual, frames can be set at 100/200/500/1000
- Display Mode: B, B/M, Color Doppler Version – B+Color, B+PDI, B+PW
- Images or Video saved as jpg, avi, & DICOM Format
- Image Frame Rate – 18 frames/second
- Software – FREE Download from App Store



5L – Linear

B/W or Color Doppler

- 7.5 or 10.0 MHz
- Penetration Depth-up to 55 mm
- 40 mm Head Length
- 128 Probe Elements



5C – Convex

B/W or Color Doppler

- 3.5 or 5.0 MHz
- 60 mm Radius Head
- Penetration Depth up to 280 mm
- 128 Probe Elements



6C – Convex

B/W or Color Doppler

- 3.5 or 5.0 MHz
- Penetration Depth-up to 305 mm
- 60 mm Radius Head
- 192 Probe Elements



Dual Headed Probe 5C – Convex & 5L – Linear

B/W or Color Doppler

Convex Probe:

- 3.5 or 5.0 MHz
- 60 mm Radius Head
- Penetration Depth-up to 280 mm
- 128 Probe Elements

Linear Probe

- 7.5 or 10.0 MHz
- 40 mm Head Length
- Penetration Depth-up to 55 mm
- 128 Probe Elements



Dual Headed Probe 5M – Micro Convex & 5L – Linear

B/W or Color Doppler

Micro Convex Probe:

- 3.5 or 5.0 MHz
- 20 mm Radius Head
- Penetration Depth-up to 280 mm
- 128 Probe Elements

Linear Probe:

- 7.5 or 10.0 MHz
- 40 mm Head Length
- Penetration Depth-up to 55 mm
- 128 Probe Elements



Dual Headed Probe w/ Phased Array 5C - Convex & 5L Linear



B/W or Color Doppler

Convex Probe

- 3.5 or 5.0 MHz
- 60 mm Radius Head
- Penetration Depth-up to 305 mm
- 128 Probe Elements

Linear Probe

- 7.5 or 10.0 MHz
- 40 mm Head Length
- Penetration Depth-up to 100 mm
- 128 Probe Elements

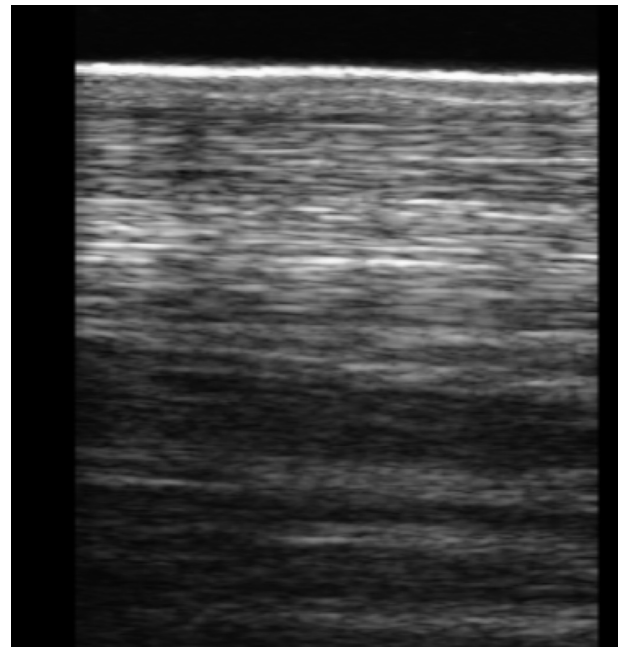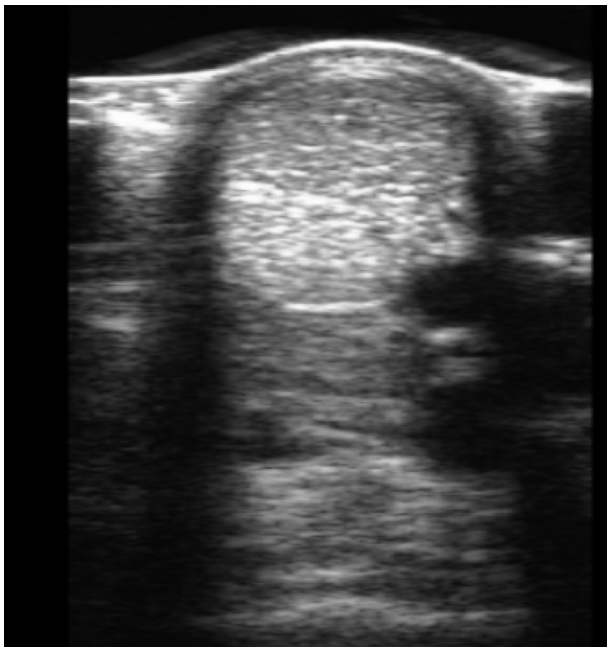

QUALITY | INTEGRITY | VALUE...privately owned since 1983

scott@VetXray.com

1.800.346.9729

www.VetXray.com