

Mixed Animal X-Ray Package

X-Ray System for the Office or Field

Digital radiography solution offering the adaptability essential for your imaging needs both in the office and in the field, complete with DR Flat Panel, portable x-ray unit, and mobile examination table.



Affordable

An affordable option for upgrading from traditional film or CR systems to digital radiography.



14" x 17" Cassette-Size

Maximized full field of view offers large format x-ray imaging for the highest quality detail in all patients without worrying about panel rotation.




Advanced Image Processing

Intelligent algorithms automatically optimize the image and provide sharp edges without noise amplification.



ULTRA Powerful

Powerful 120kV / 40mA output for small or large animal x-rays in the office or field.

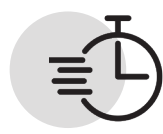

ULTRA Fast

Routine radiographs in less than .05 seconds.



ULTRA Accurate

Dual laser pointers for precise alignment, saving time, and eliminating uncertainty.



**DIAGNOSTIC
IMAGING SYSTEMS**

QUALITY | INTEGRITY | VALUE...privately owned since 1983

www.vetxray.com

Contact us:

605.341.2433

scott@vetxray.com

2325 E. Saint Charles St.
Rapid City, SD 57703





With advanced VOYANCE acquisition software and PACS compatibility, this detector has a user-friendly interface and delivers reliable performance in any environment and unexpected situation.

Conversion Screen (Scintillator)	CsI
Detector Type	TFT Direct Deposited Cesium Iodide
Pixel Pitch	140 µm
Image Data (AD conversion)	14 bit
Image Resolution	2,560 x 3,072
Image Size	14" x 17"
Data Interface	Wireless
X-Ray Sync	Lossless AED
Data Interface	1 Giga bps / IEEE 802.11n (2.4/5GHz)
Battery Performance	4hr @ 180 exposures per hr. (varies w/ standby time)
Weight	6.5 lbs

ULTRA 12040Hf

The most powerful, lightweight, high-frequency compact sized portable on the market. Designed for both small and large animals in the clinic or in the field, it produces 40 to 120 kV and up to 40 mA.

Max Output Power	2.8kW	
Input Power	Voltage	100V-240V
	Phase & Frequency	Single / 50/60 Hz
MAs Range	0.4mAs to 100mAs, 25steps	
Max. kV Deviation	±3 %	
Max. mAs Deviation	±5 %	
Display	KV / mAs: 7-segment LED	
Weight	29 lbs	

VV200

The first integrated mobile veterinary x-ray examination table to provide x-ray capabilities in examination rooms, then wheeled directly into surgery suites, with easy and accurate positioning, less movement of patient, ensuring consistent, quality results.

Table Size	56" L x 24" W top x 33" H
Table Shelf	54" L x 23" W
Weight Capacity	250 lbs
X-Ray Mounting	Easily removable portable x-ray unit
Cassette Tray	18x18 chrome brushed metal
Vertical SID	26" - 36"
X-Ray Beam Travel Movements	* Horizontal with Simultaneous Cassette Cabinet * Transverse * Oblique * Rotation * Vertical * Lateral * Angulation
Warranty	6 months



**DIAGNOSTIC
IMAGING SYSTEMS**

QUALITY | INTEGRITY | VALUE...privately owned since 1983

www.vetxray.com

Contact us:

605.341.2433
scott@vetxray.com
2325 E. Saint Charles St.
Rapid City, SD 57703

